

Nova Podiatry Workstation Cabinet - 2 Drawers

Three drawers with dedicated slide-out drill shelf

ALGEOS Order Code: NV3230B

A compact, mobile clinic cabinet designed to organise podiatry equipment and consumables while keeping a drill unit accessible at the point of care. The open equipment bay supports cable routing through the rear access point, and the three drawers provide practical storage for instruments and supplies.

The raised-edge worktop helps contain small spills and loose items, while lockable castors allow the cabinet to be positioned securely beside a treatment couch or clinical workstation.



3 drawers Including one extra-deep drawer	Slide-out drill shelf Equipment positioned within easy reach	Mobile and stable Smooth castors with locking function
Cable management Rear access aperture in equipment bay	Hygienic worktop Raised edges and easy-clean finish	Clinic-ready finish White lacquered cabinet with chrome handles

Key applications

- Podiatry and chiropody clinics
- Beauty and nail treatment environments using rotary equipment
- Medical and allied-health treatment rooms requiring compact mobile storage


Technical Specifications

Dimensions, construction and equipment accommodation

Order code	NV3230B
Product name	Nova Clinic Storage Cabinet with 3 Drawers and Drill Slider
Overall dimensions	84.5 x 50 x 37 cm (H x W x D) Inc. castor height
Cabinet body height	Approximately 78 cm before castors, based on the segmented side elevation
Open drill/equipment bay height	22.5 cm
Drill bay width and depth	The external cabinet width is 50 cm and depth is 37 cm; internal usable dimensions 47 cm and depth is 35 cm approx..
Drawers	Three drawers, including one extra-deep lower drawer
Equipment shelf	Slide-out shelf beneath the open drill bay
Cable access	Rear circular cable-management aperture
Mobility	Four smooth castors; lockable positioning feature
Worktop	Raised side and rear edges; easy-clean surface
Construction	Durable compact-laminate / lacquered cabinet construction
Colour and finish	White lacquered finish with chrome-effect handles
Products included	Cabinet only. Drill, handpiece, magnifier lamp and other accessories are not included.

Dimensions and Drill Compatibility

Use the measurements below to assess whether your drill unit can be accommodated safely.

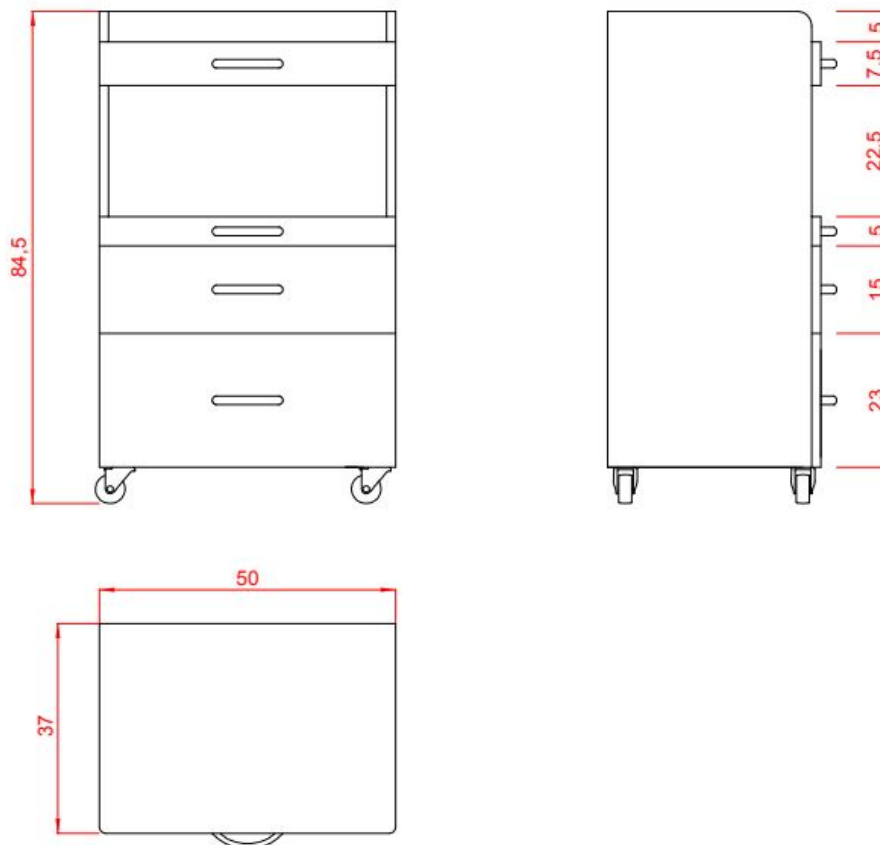

Drill handpiece holder / bracket

A dedicated drill handpiece holder is **not included**. The handpiece shown in the product photograph can be purchased separately (code SU5010)

Where a handpiece holder is used, only a compatible accessory that does not obstruct ventilation, drawer travel, castor operation or cable routing should be used. Do not drill into or permanently modify the cabinet as this may damage the finish and affect warranty and structural integrity.

Magnifier lamp attachment

A separate work-bench bracket for mounting a compatible magnifying lamp is available to purchase from Algeos, different options available depending on the type of lamp.

Product Views and Care

Illustrations show the cabinet configuration and typical equipment placement.



Cleaning and routine care

- Clean the worktop and cabinet surfaces with a soft cloth and a mild, non-abrasive cleaner suitable for lacquered furniture.
- Avoid scouring pads, solvents and aggressive chemicals that may dull or damage the finish.
- Disconnect electrical equipment before cleaning the bay or cable-access area.
- Check castors, drawer operation and handles periodically. Remove the cabinet from use if any component becomes loose or damaged.
- Maintain adequate ventilation around the drill in accordance with the drill manufacturer's instructions.

Disclaimer

Products and accessories shown in photographs are for illustration only and are not included. Specifications are based on information provided by the factory at the time of issue.