

ALGEOS INSTRUMENTS

AFFORDABLE • PROFESSIONAL • QUALITY



NEW ALGEOS INSTRUMENTS

Our newly devised range of Algeos Instruments includes everything from diabetic nippers, petite and patient-held nippers, to blacks files, casting scissors, and dressing scissors. All re-usable instruments have been forged from German Stainless Steel, have been beautifully finished, and competitively priced for the market.

Contents

Diabetic Nippers: (3) - (4)

Patient-held Nippers: (5)

Petite Nippers: (5) - (6)

Ingrown Nippers: (7) - (8)

Starter Kit: (8)

General Purpose Nippers: (9) - (11)

Cuticle, Splitter & Replacement spring: (12)

Reusable Scissors: (13) - (16)

Reusable Forceps: (17) - (19)

Blacks Files: (19) - (21)

Probes: (21)

Foot & Nail Dressers: (22)

Rasps, Splitters & Elevators: (23)

Casting Scissors: (24)

Foot Files & Applicator: (25)

Autoclave: (26)

Ordering Information: (27)

Navigation Symbols



Autoclavable



German
Stainless Steel



Safety Lock



Perfect for Retail

Instrument Care

The following chart can help you identify and prevent or remove many of the stains that typically show up on instruments during the cleaning process.

Stain/Spot Colour	Possible Causes	Solutions
Brown/Orange	Bio-burden (tissue, blood, etc.) left on nipper. Detergent pH is too high (>8). Soap/detergent residue on drapes/towels used in wrapping.	Use an eraser/buffer to rub off the stain. Choose a cleaning solution with a neutral pH. Review laundry protocol to improve soap-free rinse.
Rust	Mixing instruments made of different metals in the same cleaning or sterilising cycle. Tap water has high mineral content.	Separate instruments by metal type for cleaning and sterilisation. Use distilled water and thoroughly dry nippers after rinsing.
Black/Brown, with Pitting	Cold soaking of nippers. Detergent pH is too low (<6).	Eliminate cold soaking and potential exposure to chemicals during that process. Choose a cleaning solution with a neutral pH.
Blue/Black	Mixing instruments made of different metals in the same cleaning or sterilising cycle	Separate instruments by metal type for cleaning and sterilisation.
Blue/Grey	Improper cold sterilisation. Chemical or mineral residue.	Check solution instructions and follow recommended temperature and soak times. Use distilled water and thoroughly dry nippers after rinsing.
Spotting	Slow or improper drying of nippers. Mineral residue.	Review autoclave manufacturer's instructions. Use distilled water and thoroughly dry nippers after rinsing.
Multi/Rainbow Colour	Excessive heat during sterilising cycle.	Review autoclave instructions.

DIABETIC NIPPERS

Algeos' re-usable Diabetic Nippers have been specifically and ergonomically engineered for use on patients with diabetes. All pointed edges have been rounded off to offer a superb smooth finish, and each design is significantly lighter than most standard nippers to help prevent user fatigue and repetitive strain injury.

Diabetic General Purpose Podiatry Nipper

Diabetic General Purpose Nipper with double spring action and safety lock.

- **Joint Type:** Box.
- **Jaw:** Concave.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0045	14.5cm	£19.95



Diabetic General Purpose Podiatry Nipper

Diabetic General Purpose Nipper with double spring action and safety lock.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0046	14.5cm	£19.95



DIABETIC NIPPERS

Diabetic Fine Pointed Ingrown Podiatry Nipper

Diabetic Fine Pointed Ingrown Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0863	13cm	£18.95



Diabetic General Purpose Podiatry Nipper

Diabetic Fine Pointed Ingrown Straight Jaw Ingrown Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0864	11.5cm	£18.95



Flame Pointed Podiatry Nipper

Flame Pointed Straight Jaw Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Flame pointed.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0866	13cm	£13.95



PATIENT-HELD NIPPERS

Patient-held Nippers are designed to be single use, but they can be used twice if the same person is using them. The idea behind patient-held nippers is that the clinician uses them on the patient in clinic, and the patient is able to take them away following treatment to use again at home.

Patient-held Pedicure Nipper

Patient-held Pedicure Straight Jaw Nipper with single spring action.

- **Joint Type:** Screw.
- **Jaw:** Straight.
- **Spring:** Single spring action.
- **Handles:** Chequered.



Code	Size	Pack Size	Price
NV0870	12cm	Single	£4.95
NV0870-5	12cm	Pack of 5	£21.75

PETITE NIPPERS

Algeos has introduced reusable petite nippers into the newly devised instrument range. Petite nippers are miniature versions of the more popular full size version nipper. These would be ideally suited for women or users with smaller hands.

Petite General Purpose Podiatry Nipper

Petite General Purpose Straight Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Straight.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0868	12cm	£13.95

PETITE NIPPERS

Petite General Purpose Podiatry Nipper

Petite General Purpose Concave Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Straight.
- **Jaw:** Concave.
- **Spring:** Double spring action.
- **Handles:** Smooth.




Code	Size	Price
NV0869	12cm	£13.95



Petite Cantilever Podiatry Nipper

Petite Cantilever Concave Jaw Nipper with barrel spring action and safety lock.

- **Joint Type:** Screw.
- **Jaw:** Concave.
- **Spring:** Barrel spring action.
- **Handles:** Chequered.

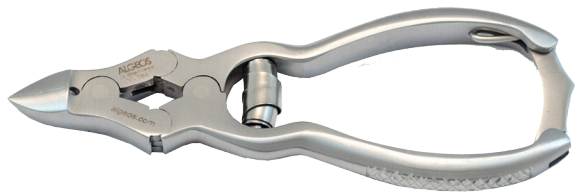

Code	Size	Price
NV0867	12cm	£22.95



INGROWN NIPPERS

Fine Pointed Ingrown Podiatry Nipper

Fine Pointed Concave Jaw Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Concave.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0011	13cm	£13.95



Fine Pointed Ingrown Podiatry Nipper

Fine Pointed Ingrown Straight Jaw Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0001	11.5cm	£13.95



Fine Pointed Ingrown Podiatry Nipper

Fine Pointed Ingrown Straight Jaw Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0002	13cm	£13.95



INGROWN NIPPERS & STARTER KIT

Fine Pointed Ingrown Podiatry Nipper

Fine Pointed Straight Jaw Nipper With double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0028	15cm	£14.95

Algeos Instrument Starter Kit

This Podiatry Instrument Starter Kit includes a selection of Algeos' most popular podiatry instruments.

This is the perfect compact kit to get you started in practice, or as a convenient instrument set to supplement your existing range. This autoclavable set is ideal for home visits.

Kit contains:

- 1 Autoclavable instrument tray and instruments with lid. (*Instrument tray measures: 195mm x 135mm x 23mm*).
- 1 14cm General Purpose Nipper (*with concave jaw*).
- 1 Black's File for the treatment of ingrown nails.
- 1 Stainless-steel Blade Handle no. 3. (*Blade handle is to be used with Nova blades only.*)
- 1 DiamonDeb Style 6" foot dresser for the reduction of hard skin and nails.
- 2 Replacement nipper springs.

Saving of **£11.80!**



Code	Price
NV0005	£19.95

GENERAL PURPOSE NIPPERS

General Purpose Podiatry Nipper

General Purpose Straight Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Chequered.



Code	Size	Price
NV0018	14.5cm	£14.95



General Purpose Podiatry Nipper

General Purpose Concave Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Box.
- **Jaw:** Concave.
- **Spring:** Double spring action.
- **Handles:** Chequered.



Code	Size	Price
NV0019	14.5cm	£14.95



General Purpose Podiatry Nipper

General Purpose Straight Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0020	14.5cm	£14.95



GENERAL PURPOSE NIPPERS

General Purpose Podiatry Nipper

General Purpose Concave Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Box.
- **Jaw:** Concave.
- **Spring:** Double spring action
- **Handles:** Smooth.



Code	Size	Price
NV0021	14.5cm	£14.95




General Purpose Podiatry Nipper

General Purpose Concave Jaw Nipper with barrel spring action and safety lock.

- **Joint Type:** Lap.
- **Jaw:** Concave.
- **Spring:** Barrel spring action.
- **Handles:** Grooved.



Code	Size	Price
NV0027	14cm	£14.95

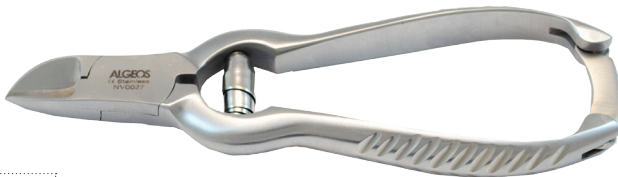


Square Handled General Purpose Podiatry Nipper

Square Handled Concave Jaw Nipper with single spring action and safety lock.

- **Joint Type:** Lap.
- **Jaw:** Concave.
- **Spring:** Single spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0017	13.5cm	£14.95

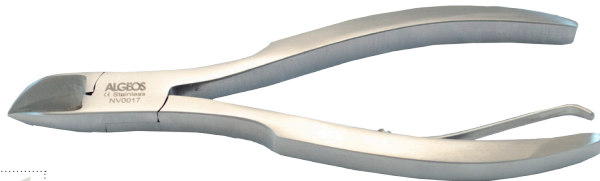

GENERAL PURPOSE NIPPERS

Cantilever General Purpose Podiatry Nipper

Cantilever Concave Jaw Nipper with double spring action and safety lock.

- **Joint Type:** Screw.
- **Jaw:** Concave.
- **Spring:** Double spring action.
- **Handles:** Chequered.



Code	Size	Price
NV0029	15cm	£24.95



Cantilever General Purpose Podiatry Nipper

Cantilever Straight Jaw Nipper with barrel spring action and safety lock.

- **Joint Type:** Screw.
- **Jaw:** Straight.
- **Spring:** Barrel spring action.
- **Handles:** Chequered.



Code	Size	Price
NV0030	15cm	£24.95



Cantilever Podiatry Nipper

General Purpose Cantilever Concave Jaw Nipper with barrel spring action and safety lock.

- **Joint Type:** Screw.
- **Jaw:** Concave.
- **Spring:** Barrel spring action.
- **Handles:** Chequered.



Code	Size	Price
NV0025	15cm	£24.95



CUTICLE NIPPER

Cuticle Nipper

Cuticle Straight Jaw Nipper with double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight (5mm).
- **Spring:** Double Spring Action.
- **Handles:** Smooth.



Code	Size	Price
NV0009	10cm	£10.95

NAIL SPLITTER

Thwaites Nail Splitter

Thwaites Nail Splitter with Straight Jaw and double spring action.

- **Joint Type:** Box.
- **Jaw:** Straight.
- **Spring:** Double spring action.
- **Handles:** Smooth.



Code	Size	Price
NV0006	14cm	£14.95

REPLACEMENT SPRING

Replacement Podiatry Nipper Springs

Algeos replacement spring - generic.



Code	Price
NV0050	£1.95



REUSABLE SCISSORS

Dressing Scissors

Dressing Scissors with blunt/sharp blades.

- **Blade Definition:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0060	16.5cm	£6.95



Dressing Scissors

Dressing Scissors with blunt/sharp blades.

- **Blade Definition:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0064	18cm	£6.95



Dressing Scissors

Dressing Scissors with blunt/blunt blades.

- **Blade Definition:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0067	13cm	£4.95



REUSABLE SCISSORS

Dressing Scissors

Dressing Scissors with blunt/sharp blades.

- **Blade Definition:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0069	14cm	£5.95

Dressing Scissors

Dressing Scissors with blunt/sharp blades.

- **Blade Definition:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0070	15cm	£5.95

Mayo Dissecting Scissors

Mayo Dissecting Scissors with straight blade.

- **Blade Definition:** Straight.
- **Handles:** Smooth/rounded.



Code	Size	Price
NV0072	14cm	£5.95

REUSABLE SCISSORS

Mayo Dissecting Scissors

Mayo Dissecting Scissors with curved blade.

- **Blade Definition:** Curved.
- **Handles:** Smooth/rounded.



Code	Size	Price
NV0073	14cm	£5.95



Metzenbaum Scissors

Metzenbaum Scissors with straight blade.

- **Blade Definition:** Straight.
- **Handles:** Smooth/rounded.



Code	Size	Price
NV0076	14cm	£5.95



Lister Bandage Scissors

Lister Bandage Scissors with straight angled jaw.

- **Blade Definition:** Straight/angled.
- **Handles:** Smooth/rounded.



Code	Size	Price
NV0078	14cm	£5.95



REUSABLE SCISSORS & FORCEPS

Utility Scissors

Utility Scissors with serrated blade and plastic handles.

- **Blade Definition:** Straight/serrated blade.
- **Handles:** Smooth/flat.



Code	Size	Colour	Price
NV0086	18cm	Black	£2.95
NV0135	18cm	Blue	£2.95

Halstead Mosquito Forcep

Halstead Mosquito Forcep with serrated jaw.

- **Joint Type:** Box.
- **Jaw:** Serrated.
- **Blades:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0089	12.5cm	£3.95

Spencer Wells Artery Forcep

Spencer Wells Artery Forcep with serrated jaw.

- **Joint Type:** Box.
- **Jaw:** Serrated.
- **Blades:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0091	12.5cm	£3.95

REUSABLE FORCEPS

Adson Dissecting Forcep

Adson Dissecting Forcep with straight jaw.

- **Jaw:** Toothed/straight.
- **Handles:** Smooth/flat with grooved body.



Code	Size	Price
NV0098	12cm	£3.95

Adson Dissecting Forcep with 1:2 Teeth

Adson Dissecting Forcep with straight jaw and 1:2 teeth.

- **Jaw:** Straight.
- **Handles:** Smooth/flat with grooved body.
- **Teeth:** 1:2.



Code	Size	Price
NV0113	12cm	£2.95

Iris Dissecting Forcep

Iris Dissecting Forcep with serrated jaw.

- **Jaw:** Serrated.
- **Handles:** Smooth/flat with grooved body.
- **Pin Mechanism:** Yes.



Code	Size	Price
NV0104	13cm	£2.95

REUSABLE FORCEPS

Rampléy Sponge Forcep

Rampléy Sponge Forcep with serrated jaw.

- **Joint Type:** Box.
- **Jaw:** Serrated.
- **Blades:** Straight.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0111	24cm	£6.95




Lister Sinus Forcep

Lister Sinus Forcep with serrated jaw.

- **Joint Type:** Box.
- **Jaw:** Serrated.
- **Handles:** Smooth/flat.
- **Blades:** Straight.
- **Non Locking.**



Code	Size	Price
NV0100	15cm	£4.95

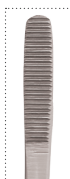

Kocher Forcep

Kocher Forcep with serrated jaw.

- **Joint Type:** Box.
- **Jaw:** Serrated.
- **Handles:** Smooth/flat.
- **Teeth:** 1:2.



Code	Size	Price
NV0118	13cm	£4.95



FORCEPS & BLACKS FILES

Blade Removing Forcep

Blade Removing Forcep with serrated jaw and curved blade.

- **Joint Type:** Lap.
- **Jaw:** Serrated.
- **Blades:** Curved.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0120	16cm	£4.95



Blacks File

The stainless steel fine Blacks file is designed for the treatment of ingrown nails.

- **Shape/End:** Straight/single-ended.
- **Handles:** Patterned hexagonal.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0126	13cm	£2.95



Blacks File

The stainless steel standard/medium Blacks file is designed for the treatment of ingrown nails.

- **Shape/End:** Straight/single-ended.
- **Handles:** Patterned hexagonal.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0128	13cm	£2.95



BLACKS FILES

Blacks File

The stainless steel Blacks file micro head is designed for the treatment of ingrown nails.

- **Shape/End:** Straight/single-ended.
- **Handles:** Patterned hexagonal.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0129	13cm	£2.95

Blacks File

The stainless steel Blacks file super cut head is designed for the treatment of ingrown nails.

- **Shape/End:** Straight/single-ended.
- **Handles:** Patterned hexagonal.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0130	13cm	£2.95

Blacks File

The stainless steel, double ended fine pointed Blacks file is designed for the treatment of ingrown nails.

- **Shape/End:** Straight/double-ended/fine pointed.
- **Handles:** Patterned hexagonal.
- Rigid grip to increase comfort and prevent slippage.




Code	Size	Price
NV0131	14cm	£3.95

BLACKS FILES & PROBES

Blacks File

The stainless steel, double ended straight angled Blacks file is designed for the treatment of ingrown nails.

- **Shape/End:** Angled/straight/double-ended with medium cut.
- **Handles:** Patterned hexagonal.
- Rigid grip to increase comfort and prevent slippage.

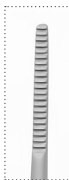

Code	Size	Price
NV0132	14cm	£3.95

Probe

Probe with flat end.

- **Shape/End:** Straight/flat/single-ended.
- **Handles:** Round/ribboned.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0138	13cm	£2.95

Probe

Probe with spoon end.

- **Shape/End:** Straight/spoon/single-ended.
- **Handles:** Round/ribboned.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0140	13cm	£2.95

FOOT & NAIL DRESSERS

DiamonDeb Style Foot Dresser

Professional DiamonDeb Style foot dresser for the filing of nails and hard skin.

- **Autoclavable:** 20cm - Yes.
15cm - No.
- Bulk packs available.



Code	Size	Price
NV0143	15cm	£1.95
NV0144	20cm	£3.25

Packs of 5

Code	Size	Price
NV0143-5	15cm	£7.00
NV0144-5	20cm	£14.95



DiamonDeb Style Foot Dresser

Professional DiamonDeb Style foot dresser for the filing of nails and hard skin.



Code	Size	Price
NV0871	15cm	£3.25



DiamonDeb Style Nail Dresser

Professional DiamonDeb Style nail dresser shapes nails smoothly and helps prevent splitting and peeling.



Code	Size	Price
NV0145	20cm	£1.95
NV0146	15cm	£1.95



NAIL RASP, SPLITTER & ELEVATOR


Nail Rasp

Double-sided Nail Rasp with flat handle.

- **Handles:** Smooth/flat.




Code	Size	Price
NV0147	17cm	£3.95

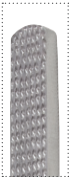

Francis Nail Splitter

Francis Nail Splitter with flat patterned handle.

- **Handles:** Flat/patterned/grooved body.

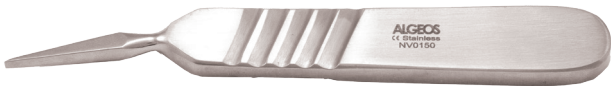

Code	Size	Price
NV0150	11.5cm	£3.95



Nail Locke Elevator

Nail Locke Elevator with round knurled handle.

- **Shape/End:** Straight/single-ended.
- **Handles:** Round/knurled.
- Rigid grip to increase comfort and prevent slippage.



Code	Size	Price
NV0153	13cm	£2.95



CASTING SCISSORS

Wire Cutting Scissors

Stainless Steel Wire Cutting Scissor with serrated blades.

- **Jaw:** Angled.
- **Blade:** Serrated.
- **Handles:** Smooth/flat.



Code	Size	Price
NV1126	12.5cm	£4.95

Wire Cutting Scissors

Stainless Steel Wire Cutting Scissor with tungsten carbide blades.

- **Jaw:** Straight.
- **Blade:** Tungsten Carbide.
- **Handles:** Gold plated/flat/smooth.



Code	Size	Price
NV0085	12.5cm	£11.95

Lister Bandage Scissors

Stainless Steel Lister Bandage Scissor with tungsten carbide blades.

- **Jaw:** Angled.
- **Blade:** Tungsten carbide.
- **Handles:** Smooth/flat.



Code	Size	Price
NV0075	18cm	£14.95

FOOT FILES & APPLICATOR

Plastic Foot File

Disposable Plastic Foot File that helps remove dry, dead skin from heels.

- Ensures a soft, smooth skin finish.
- Features a large file area.
- Curved handle for easy use.
- Double sided.



Code	Description	Price
NV0148	Plastic Foot File (Singles)	£0.95

Disposable Foot Files

A cost effective alternative to reusable foot files.

- Non-sterile, single patient use files.
- Double sided.
- **Pack Size:** Pack of 10.



Code	Price
NV0155	£2.50

Tubular Gauze Applicator

The shape of a sprung sugar tong for quick and easy application of tubular gauze or tubular bandages.

- **Size:** 9cm.



Code	Description	Price
NV0161	Small	£2.95
NV0162	Large	£2.95

AUTOCLAVE

MELAtronic® 15 EN+ 'Class S' (7L) Autoclave

Product Features & Benefits:

- Compact dimensions - less installation space required.
- User-friendly operation.
- Quick operation times.
- Effective drying.
- Compatible with MELAflash CF-Card or MELAprint 42 Log Printer for perfect documentation (separate product information available).
- Integrated water storage tank* - no water connection required.
- Conforms to European Standard EN13060.
- Four fully automatic, microprocessor controlled programme processes.
- Internal space for three trays (sold separately)*.
- Accurate presentation of the parameters on the digital display.
- Short sterilisation time with additional drying or preheat functions.
- Automatic water supply and water quality control.
- Electric door lock with mechanical interlocking.
- Tray mount (holds 3 trays) included as standard*.
- Made in Germany.

*Distilled water (from a water distiller) is required.



Code	Description	Price
ML10615	MELAtronic 15EN+ Autoclave + Professional Installation	£2,740.00
ML00150	MELAG Autoclave Tray (35cm x 12cm) for MELAtronic 15EN+	£36.95
ML70670	Foil Holder for MELAtronic 15EN+ (for upright sterilisation of sealed instruments - for better drying results)	£56.95

Autoclavable Trays

- Autoclavable up to 135°C.
- Can be sterilised by autoclave, ethylene oxide or hydrogen plasma.
- Manufactured from high-quality polypropylene.
- Improved removal of blood residue compared to stainless steel.



Code	Size/Type	Price
NV0822	L195mm x W135mm x H23mm	£3.25
NV0823	L270mm x W150mm x H30mm	£3.25
NV0821	Lid for NV0823	£3.25
NV0824	Lid for NV0822	£3.25

Visit www.algeos.com
for the full range of
autoclaves and
accessories.

ORDERING INFORMATION

Our friendly, knowledgeable, and experienced Customer Service team is available Monday to Friday from 8:30am to 5pm.

 **Telephone:** +44 (0)151 448 1228

 **E-mail:** orders@algeos.com

 **Fax:** +44 (0)151 448 1008

 **Post:** A. Algeo Ltd
Sheridan House
Bridge Industrial Estate
Speke Hall Road
Liverpool, L24 9HB

 **Online:** www.algeos.com

Our complete online catalogue is available **24-7** enabling you to browse at your leisure and place orders when it's most convenient.





Algeos Head Office

Sheridan House, Bridge Industrial Estate,
Speke Hall Road, Liverpool, L24 9HB

t: +44 (0) 151 448 1228

f: +44 (0) 151 448 1008

e: info@algeos.com

w: algeos.com